

## **Steriliser**


# SN75plus


Safer than disinfection, hot-air sterilisation of containers and instruments



The indispensable safety feature for this hot air steriliser: Setpoint Wait. This means that the programme for sterilisation only starts when the set temperature has been reached. This feature can also be used with freely positionable Pt100 temperature sensors. Here the sterilisation time only begins when the set temperature has been reached at all measurement points, and reliable sterilisation is guaranteed at all times.

On this page, you can find all the essential technical data on the Memmert hot air steriliser. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at <a href="mailto:sales@memmert.com">sales@memmert.com</a>.

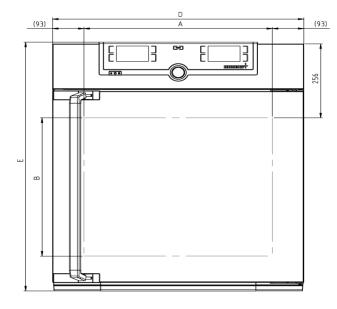

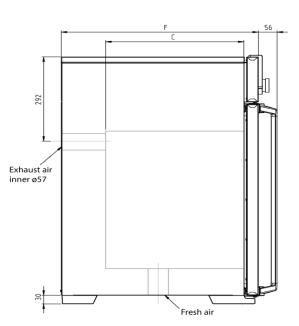

Temperature		
Setting temperature range	+20 to +250 °C	
Setting accuracy temperature	up to 99.9 °C: 0.1 / from 100 °C: 0.5	
Working temperature range	at least 5 above ambient temperature to +250 °C	
Temperature sensor	2 Pt100 sensors DIN Class A in 4-wire-circuit for mutual monitoring, taking over functions in case of an error	
Control technology		
ControlCOCKPIT	TwinDISPLAY. Adaptive multifunctional digital PID-microprocessor controller with 2 high-definition TFT-colour displays.	
Language setting	German, English, Spanish, French, Polish, Czech, Hungarian	
Timer	Digital backwards counter with target time setting, adjustable from 1 minute to 99 days	
Function HeatBALANCE	adapting the distribution of the heating performance of the upper and lower heating circuit from -50 $\%$ to +50 $\%$	
Function SetpointWAIT	the process time does not start until the set temperature is reached	
Calibration	three freely selectable temperature values	
adjustable parameters	temperature (Celsius or Fahrenheit), air flap position, programme time, time zones, summertime/wintertime	
Ventilation Convection Fresh air admixture	natural convection  adjustment of pre-heated fresh air admixture by air flap control in 10 % steps for each segment	
Vent	individually  vent connection with restrictor flap	
Communication  Documentation	programme stored in case of power failure	
Programming	AtmoCONTROL software on a USB stick for programming, managing and transferring programmes via Ethernet interface or USB port	
Safety		
Temperature control		
Temperature control	mechanical temperature limiter TB, protection class 1 according to DIN 12880 to switch off the heating approx. 20°C above nominal temperature	
Temperature control		
·	approx. 20°C above nominal temperature  overtemperature monitor TWW, protection class 3.1 or adjustable temperature limiter TWB, protection	
Temperature control	approx. 20°C above nominal temperature  overtemperature monitor TWW, protection class 3.1 or adjustable temperature limiter TWB, protection class 2, selectable on display  additionally integrated over- and undertemperature monitor "ASF", automatically following the setpoint value at a preset tolerance range, alarm in case of over- or undertemperature, heating is switched off	

•			
Stan	dard	eauip	mant
Juan	uaiu	Cuuin	,,,,,

Door	fully insulated stainless steel door with 2-point locking (compression door lock)
Internals	2 stainless steel grid(s), electropolished
Works calibration certificate	Calibration at +160°C

## Stainless steel interior

Interior	easy-to-clean interior,made of stainless steel, reinforced by deep drawn ribbing with integrated and protected large-area heating on four sides
Volume	74
Dimensions	w <sub>(A)</sub> x h <sub>(B)</sub> x d <sub>(C)</sub> : 400 x 560 x 330 mm
Max. number of internals	6
Max. loading of chamber	120 kg
Max. loading per internal	20 kg

# Textured stainless steel casing

Dimensions	w <sub>(D)</sub> x h <sub>(E)</sub> x d <sub>(F)</sub> : 585 x 944 x 514 mm (d +2x56mm door handle)
Housing	rear zinc-plated steel

#### **Electrical data**

Voltage	230 V, 50/60 Hz
Electrical load	approx. 2500 W
Voltage	115 V, 50/60 Hz
Electrical load	approx. 1800 W

## **Ambient conditions**

Set Up	The distance between the wall and the rear of the appliance must be at least 15 cm. The clearance from the ceiling must not be less than 20 cm and the side clearance from walls or nearby appliances must not be less than 5 cm.
Altitude of installation	max. 2,000 m above sea level
Ambient temperature	+5 °C to +40 °C
Humidity rh	max. 80 %, non-condensing
Overvoltage category	II
Pollution degree	2

# Packing/shipping data

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-RegNo.	DE 66812464
Dimensions approx incl. carton	w x h x d: 730 x 1130 x 670 mm
Net weight	approx. 66 kg
Gross weight carton	approx. 85 kg

# Standard units are safety-approved and bear the test marks







